

Sere-med

Life Enhancing Innovations

No anesthesia needed.

Powerful. Precise.
For better treatments
and results.

—

www.seremed.com



WHAT TO TREAT

Smaller lesions on the face, nose, eyelids, ears, legs and feet are often considered not important enough for surgery under general or local anaesthetic.

CryoPen offers an innovative solution in these cases.

Pets can be in and out in just minutes and they emerge bright and bouncy, rather than having to be monitored for two to three hours.

There's no cutting or bleeding.

Therefore, nothing for your patient to lick or pull on. Minimizing risk of infection.



Tumor of the eyelid

Treatment time: 5-15 seconds



Common wart

Treatment time: 15 seconds



Tumor of the eyelid

Treatment time: 10-20 seconds



Skin tumor

Treatment time: 15-30 seconds

SOME EXAMPLES

Most common pathologies to treat are: skin tumours, eyelid tumours, warts, papilloma, cysts, perianal and sebaceous adenoma and more...



Before



After




HOW TO TREAT

Tissue freezing is quick. Treatments only take seconds, with no anaesthesia or pre- or post-op care required.

The CryoPen works with very practical disposable cartridges and emits a fine jet of nitrous oxide under high pressure.


Thanks to the device's extreme precision there is no collateral damage to healthy tissue. Load the device, choosing from the selection of applicators for treatment from a pinpoint, to 1mm, up to a width of 20mm. Point, click and freeze...

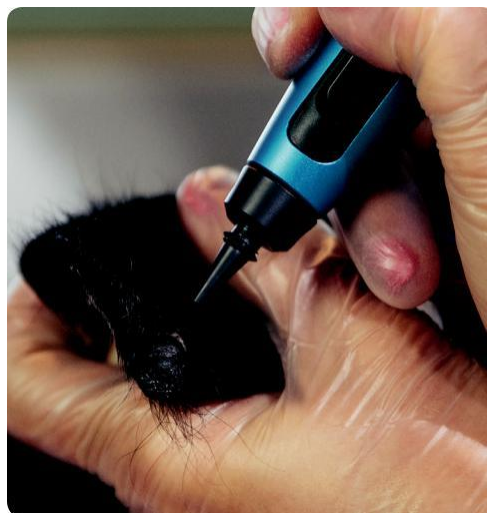

Cryosurgery, of course, is nothing new. It's recognized as THE reference treatment in dermatology. However, veterinarians who have struggled with cryoliquid spills, imprecise cotton buds and laborious deliveries and handling of hazardous liquid nitrogen will find that the CryoPen is a totally different experience.



Before

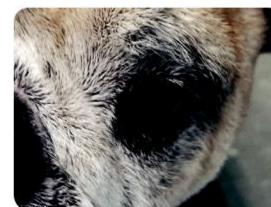

After

CryoPen B+

To be used with 23.5g cartridges.

3 applicators: Orange Ø 1-4mm - White Ø 2-6mm - Green Ø 4-8mm



CryoPen O

To be used with 8g & 16g cartridges.

3 applicators: Blue Ø 1-3mm - White Ø 2-6mm - Green Ø 4-8mm



CryoPen X+

To be used with 8g & 16g cartridges.

4 standard applicators: Blue Ø 1-3mm - White Ø 2-6mm - Green Ø 4-8mm - Yellow Ø 7-20mm
2 longer 60mm applicators: Blue Ø 1-3mm - White Ø 2-6mm



CryoPen XP

To be used with 8g & 16g cartridges.

4 applicators: Blue Ø 1-3mm - Silver Ø 2-4mm - Green Ø 3-5mm - Yellow Ø 5-8mm



Cartridges

8g - 100 seconds of constant gas flow
16g - 200 seconds of constant gas flow
23.5g - 300 seconds of constant gas flow